

# DCI II MANUAL CONTROL DENTAL DELIVERY SYSTEM FOR TWO HANDPIECES

## Instructions for Use



4102 / 4106

### GENERAL INFORMATION

#### APPLICABLE MODEL

4XXX

#### INTENDED APPLICATION AND USE

The dental delivery system is a device intended to support the instruments used by the dental practitioner, delivering those instruments to an accessible position during a dental procedure. This device may control and be the means of delivering compressed air, water, vacuum and low voltage electricity to a variety of instruments commonly used in dental practice.

### OPERATION AND FEATURES

#### CONTROLS

The *handpiece selector toggle* controls which handpiece is active. The toggle position (left or right) indicates the active handpiece. It is located on the underside of the control dental delivery system, in the center.


The *water coolant on/off toggle* turns the air signal on or off to the water coolant flow control/relay valve for handpiece #2. It is located on the underside of the control dental delivery system, in the left front.

The *water coolant flow control/relay knob* controls water to handpiece #2. It is located on the underside of the control dental delivery system, in the right front.

If a wet handpiece option is installed for handpiece #1, the flow control knob for it is installed behind the knob for handpiece #2.

The *drive air pressure gauge* indicates the operating pressure of the selected handpiece. It is located on the front top center of the control dental delivery system cover.

The *drive air adjustment screws* are located on the underside of the control dental delivery system behind the handpiece selector toggle. Handpiece drive air pressures should be adjusted to the handpiece manufacturer's recommendation. See the section on *adjustments*.

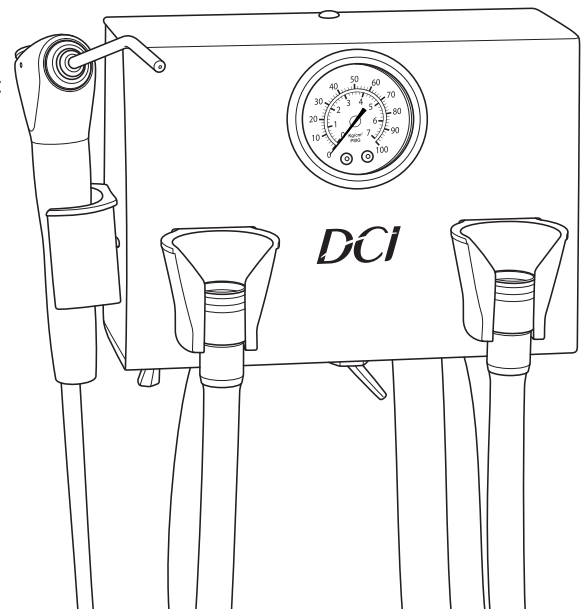

#### Cleaning and Maintenance

*Do not use powdered cleansers, scouring pads, or abrasive scrubbers on any of the finished metal surfaces in this control dental delivery system, i.e. the syringe or the foot control disc. Sodium Hypochloride will also damage these surfaces.*

#### SYRINGE

The dental delivery systems comes with a Precision Comfort syringe. The syringe is packaged in a shipping envelope, with operating instructions and a repair kit attached. The syringe goes in the holder on the far-left side of the dental delivery system.

## FOOT CONTROL

The foot control is an economy style. Handpiece speed is controlled with the foot control disc. Varying pressure on the foot control disc controls speed. Air coolant is also provided when you step on the foot control disc.

## CONTROL HEAD

The control head can be cleaned with most commonly available surface disinfectants. Do not use any Sodium Hypochloride solutions, or any cleansers containing alcohol. These may cause paint and finish discoloration.

## DENTAL DELIVERY SYSTEM WATER LINE MAINTENANCE

The Center for Disease Control and the American Dental Association can provide recommendations on when to flush your system, for how long, and what to use.

## ADJUSTMENTS



**NOTE:** All of the following adjustments should be made with a bur in the handpiece. Running a handpiece without a bur can damage the handpiece.

## HANDPIECE DRIVE AIR PRESSURE

Refer to the manufacturer's literature to determine the recommended drive air operating pressure for your handpieces.

You will need a 3/32" hex key to make these adjustments. Install a bur in the handpiece to be tested.

With the handpiece selector toggle select the handpiece you want to adjust. Place the water coolant on/off toggle in the off position. Step on the foot control disc until the handpiece is running at maximum speed.

Handpiece pressure adjustment screws are located in the rear center, behind the handpiece selector valve. Using the 3/32" ball driver adjust handpiece pressures. The front adjusting screw is for handpiece #2 and the rear adjusting screw is for handpiece #1. Turn the adjustment screw counterclockwise until the pressure gauge reads a little more than the recommended operating pressure. Then turn the screw until the pressure gauge indicates the recommended operating pressure. Repeat this step for the remaining handpiece.

## WATER COOLANT FLOW CONTROL

Place the water coolant on/off toggle in the on position. Install a bur in the handpiece to be tested. Press on the foot control disc until the handpiece is running at half operating speed. Adjust the water coolant flow control knob for handpiece #2 until a fine spray is present around the bur. Very little water coolant is required to attain the appropriate spray. Repeat process for handpiece #1 if a wet handpiece option is installed.

